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Revised in October, 2016

1. Introduction

This Hemoglobin A1c Control Set has been designed exclusively for the quality control of Tosoh Automated Glycohemoglobin Analyzers, Standard Analysis Mode and Variant Analysis Mode. This product is intended for evaluation and monitoring of the performance of HbA_{1c} assays on those analyzers by being assayed as samples.

This product consists of two levels of HbA_{1c} (% or mmol/mol) to cover the clinically significant range for diabetes control.

2. Prior to Use

Inspect the packing and the exterior of the vial for any signs of damage prior to use. If any damage is visible, contact your local Tosoh sales representative.

Confirm that the following document is included in the package.

Instructions For Use (this document) 1 copy

3. Warnings and Precautions

- 1) This product is for *in vitro* diagnostic use only.
- 2) This product is intended for use on Tosoh Automated Glycohemoglobin Analyzers (HLC-723[®]G7, HLC-723G8, HLC-723GX and HLC-723G11), Standard Analysis Mode and Variant Analysis Mode.
- 3) After reconstitution, the product must be tightly sealed and stored in an upright position at 2 - 8 °C.
- 4) The reconstituted product can be stored under freezing condition. In such a case, it must be frozen immediately after reconstitution.
- 5) Repeated freeze-thaw cycles must be avoided.
- 6) In case an erroneous result is obtained with the reconstituted and frozen product, use another vial of the product freshly reconstituted.
- 7) Human blood used in the preparation of this product has been tested by FDA-approved methods and found negative for the presence of HBsAg and antibodies to HCV and HIV-1. Because no test method can offer complete assurance that products derived from human blood will not transmit infectious agents, it is recommended that this product be handled with the same precautions as used for patient specimens.
- 8) Do not use this product beyond the expiration date.
- 9) In case of accidental ingestion, rinse the mouth and throat with excess water and immediately call for medical attention.
- 10) For safe waste disposal, it is recommended that each laboratory complies with established laboratory procedures and local, state, and federal regulations.

4. Content

Catalogue No.	Description	Package content
0021974	Hemoglobin A1c Control Set	HbA _{1c} Control Level 1 : 4 vials × 0.5 mL HbA _{1c} Control Level 2 : 4 vials × 0.5 mL

This product has been prepared from human blood cells with two significant level of HbA_{1c} (% or mmol/mol) and lyophilized. The total hemoglobin concentration is about 40 g / L (4 g / dL) after reconstitution.

5. Related Components

Hemoglobin A1c Calibrator Set Catalogue No. 0018767

Traceable to the IFCC reference method. The NGSP values are assigned using the Master equation (See Section 7 below).

6. Storage and Stability

- 1) Before opening, Hemoglobin A1c Control Set should be stored at 2 - 8 °C.
- 2) The product is stable for 7 days after opening or reconstitution provided that the vial is kept tightly sealed and refrigerated at 2 - 8 °C.
- 3) When the product is frozen after reconstitution, it should be stored at - 20 °C or below. Then, it will be stable for up to 30 days.

4) The expiration date for this product is shown below.

Expiration date : 2021-11 (YYYY-MM)

7. Assigned Values

Lot No. AB9040

	Level 1	Level 2
IFCC aligned value (mmol/mol)	30 ± 3	88 ± 5
NGSP aligned value (%)	4.9 ± 0.3	10.2 ± 0.5

IFCC : International Federation of Clinical Chemistry and Laboratory Medicine
The IFCC aligned values are traceable to the IFCC reference method.

NGSP : National Glycohemoglobin Standardization Program
The NGSP aligned values were calculated using the following conversion equation (Master equation) from the IFCC aligned values:

$$\text{NGSP (\%)} = 0.09148 \times \text{IFCC (mmol/mol)} + 2.152$$

Ref. : Geistanger A. et al. Statistical Methods for Monitoring the Relationship between the IFCC Reference Measurement Procedure for Hemoglobin A1c and the Designated Comparison Methods in the United States, Japan, and Sweden. Clin Chem 2008; 54: 1379-1385.

The assigned values are specific to each lot of the product. The observed values may vary during the life time of the product and may depend on the used instrument and reagents. Each laboratory should establish its own criteria for acceptable range of variations.

8. Preparation

- 1) Tear off the metal seal with plastic flip top and carefully remove the rubber cap.
(NOTE) A sudden rush of air can result in a loss of lyophilized material.
- 2) Using volumetric pipettes, accurately reconstitute the lyophilized material with 0.5 mL of purified water.
- 3) Place the removed cap at step 1) on the vial and let the material stand at room temperature for 30 minutes.
- 4) Swirl gently but thoroughly before use to ensure homogeneity.
- 5) Dilute the reconstituted product with purified water by 51 times (10 µL of the reconstituted product to 0.5 mL of purified water) before performing the assays.

9. Long Term Storage Procedure

When the reconstituted product is stored at - 20 °C or below, it will remain stable for up to 30 days. The reconstituted product must be frozen immediately after reconstitution.

- 1) It is recommended that the reconstituted product is frozen by being dispensed into plastic tubes.
- 2) After dispensing, the tubes should be tightly sealed and immediately frozen at - 20 °C or below.
- 3) Prior to assay, slowly bring the frozen product to room temperature, and wait until it completely thaws. After that, dilute it with purified water by 51 times (10 µL of the thawed product to 0.5 mL of purified water) and gently mix it.
- 4) Repeated freeze-thaw cycles must be avoided.

10. Assay Procedure

After reconstitution and dilution, the Hemoglobin A1c Control Set should be treated in the same manner as an unknown patient specimen manually diluted. It should be analyzed on a Tosoh Automated Glycohemoglobin Analyzer following the instructions in the Operator's Manual.

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1. Geistanger A. et al. Statistical Methods for Monitoring the Relationship between the IFCC Reference Measuring Procedure for Hemoglobin A1c and the Designated Comparison Methods in the United States, Japan, and Sweden. Clin Chem. 54: 1379 - 1385, 2008

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注册人/生产企业名称：Tosoh Corporation 东曹株式会社
住所：4560, Kaisei-cho, Shunan, Yamaguchi 746-8501, Japan
生产地址：4560, Kaisei-cho, Shunan, Yamaguchi 746-8501, Japan
联系方式：+81 3 5427 5181

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European Conformity

Manufacturer

Authorized representative in the European Community

Catalogue number / Part number

In vitro diagnostic medical device

Consult instructions for use

Use-by date / Use by

Batch code / Lot number


Temperature limitation

Net volume (after reconstitution for lyophilized material)

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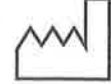This manual may not be reprinted or copied in whole or in part without written consent of Tosoh Corporation. The contents of the manual are subject to change without notice.

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"HLC" and "HLC-723" are the registered trademarks of Tosoh Corporation in Japan, etc. "G11" is the registered trademark of Tosoh Corporation.

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